


CovaWound™ Foam Border

Non-Adherent Foam Dressing with Adhesive Border

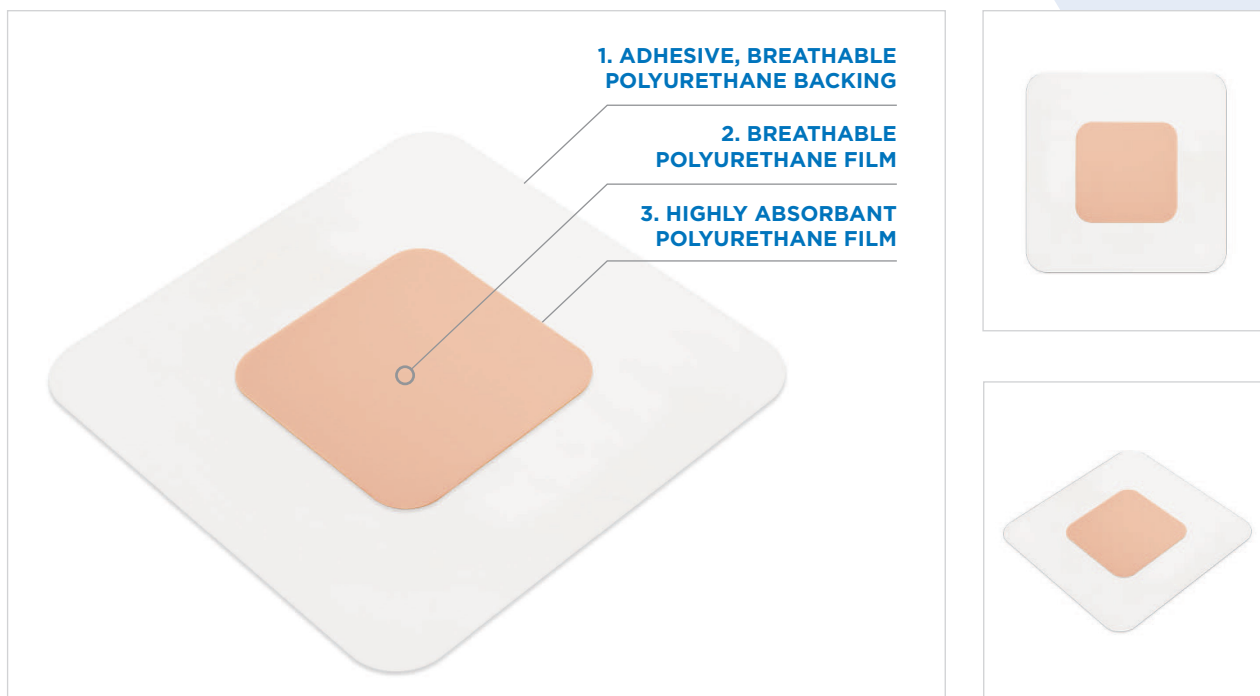

A non-adherent foam dressing with a two layer construction composed of highly absorbent polyurethane foam with a breathable polyurethane film backing. The absorbent polyurethane foam core is designed to absorb and retain exudate, while the breathable outer layer allows moisture vapour transmission away from the dressing. CovaWound™ Foam dressing enables healing to progress by creating and maintaining a moist healing environment. The breathable pink outer layer protects against eschar formation and is also waterproof, protecting against external contamination.

CovaWound™ Foam dressing is soft and conformable providing an optimal moist healing environment and excellent patient comfort. As CovaWound™ Foam dressing absorbs exudate, it may swell slightly but will not leave debris in the wound.

1. ADHESIVE, BREATHABLE POLYURETHANE BACKING

Waterproof, protecting against external contamination.

2. BREATHABLE POLYURETHANE FILM

Allows moisture vapour transmission away from the dressing.

3. HIGHLY ABSORBENT POLYURETHANE FOAM

Designed to absorb and retain exudate, while still maintaining a moist wound environment.

Exudate level

👉 👉 👉 **Moderate**

Exudate level indicates the highest level of exudate that the dressing can accommodate.

Product Codes and Sizes: Available to order through NHS Supply Chain Advanced Wound Care Framework and Alliance Healthcare

Dressing Size (cm)	Pad Size	Pack Size	Covalon Article Number	PIP Code	Catalogue Code
7.5 x 15	4 x 11	10	2030411		ELA938
7.5 x 20	4 x 16	10	2030416		ELA935
10 x 10	5 x 5	10	2010505	406-5728	ELA928
10 x 20	5 x 15	10	2010515	406-5793	ELA929
15 x 15	10 x 10	10	2011010	406-5827	ELA932
12.5 x 22.5	9 x 19	10	2010919		ELA930
15 x 15	11 x 10	10	2011110		ELA931
15 x 20	10 x 15	10	2011015	406-5744	ELA933
17 x 20	13 x 16	10	2011316		ELA934

Indications for Use

CovaWound™ Foam dressing does not adhere to the wound. It is indicated for the management of light to moderately exuding wounds such as pressure, leg, and diabetic foot ulcers, superficial and partial-thickness burns, donor sites, trauma wounds and post-operative surgical wounds.

Application

- Clean the wound area with normal saline or distilled water.
- Gently dry the surrounding skin and ensure that it is free of any residue.
- Choose the size and shape of the dressing which is most suitable, the dressing should extend 2 cm (3/4") beyond the wound margin.
- Remove the dressing from the packaging.
- Remove the release liner from the dressing and apply the white side of the non-adherent foam dressing to the wound.
- Firmly ax the adhesive border to the surrounding skin.
- Remove the back release liner and smooth out the dressing to ensure it is firmly in place.

CovaWound™ Advanced Wound Care Range



CovaWound™
Silicone Border &
Silicone Border Lite




CovaWound™
Silicone & Silicone Lite



CovaWound™
Sacral



CovaWound™
Alginate

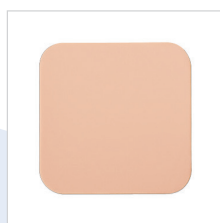

CovaWound™
Foam



CovaWound™
Adhesive Foam



CovaWound™
Hydrocolloid &
Hydrocolloid Thin/
Transparent



CovaWound™
SuperAbsorbent